

# Operating Manual for H2601 Rotary Brush Beauty Instrument







## ●Working Principle of Rotary Brush:


Using the energy to push the brush twirl for cleaning and massage skin. The function of cleaning can be wiped dirt in the surface skin, oil in the face, leftover of cosmetics, died horniness off. While the function of massage can be make skin smooth and tender.

## ●Specification, Model & Technical Data:

|                           |                         |
|---------------------------|-------------------------|
| Product                   | Beauty Instrument       |
| Model No                  | H2601                   |
| Functions                 | Rotary Brush            |
| Nominal Working Voltage   | ac. 100-130V / 220-240V |
| Nominal Working Frequency | 50/60Hz                 |
| Output Power              | 20 W                    |
| Dimension (L*W*H)         | 442*347*147mm           |
| Net Weight                | 1.8Kg                   |

## ●Packing List :
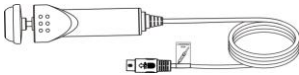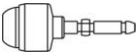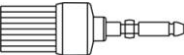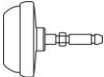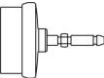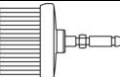
| Name                     | Graphic                                                                             | Quantity |
|--------------------------|-------------------------------------------------------------------------------------|----------|
| Host Machine             |  | 1 Round  |
| Power cable              |  | 1 Piece  |
| Brush handle             |  | 1 Piece  |
| Crystallized small brush |  | 1 Piece  |
| Fuzz small brush         |  | 1 Piece  |
| Crystallized big brush   |  | 1 Piece  |
| EVA big brush            |  | 1 Piece  |
| Fuzz big brush           |  | 1 Piece  |

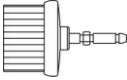
|                    |                                                                                   |         |
|--------------------|-----------------------------------------------------------------------------------|---------|
| shag big brush     |  | 1 Piece |
| Instruction Manual |                                                                                   | 1 Copy  |

## ●Host Machine:



## ●Accessories :


### Brush handle

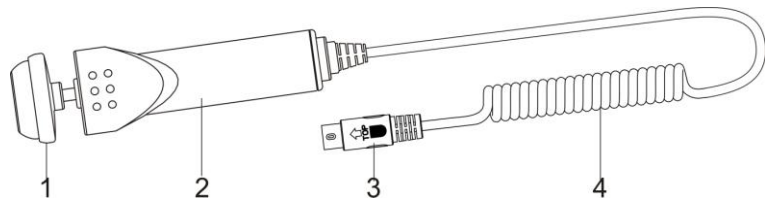


1. Brush head      2. Brush handle      3. Brush flexible wire plug      4. Brush flexible wire

### Brushes

1. Crystallized small brush: Clean small area horny and callus of hands and feet.
2. Fuzz small brush: Clean and massage small area skin of face etc..
3. Crystallized big brush: Clean large area horny and callus of hands and feet.
4. EVA big brush: Clean and massage large area skin of face etc.
5. Fuzz big brush: Clean and massage large area skin of face etc.
6. Shag big brush: Clean and massage large area skin of face and other position.


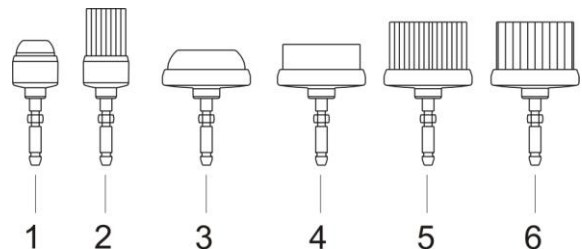


## ●Introduction to Control Panel:

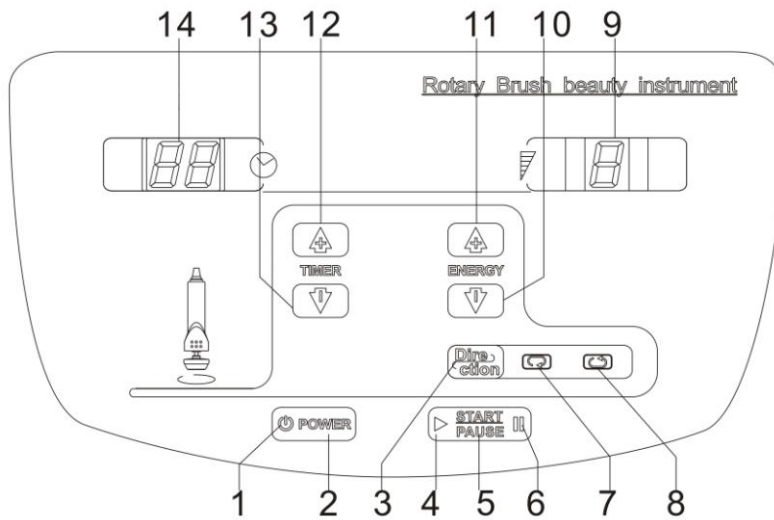

- |                                                       |                         |                                      |
|-------------------------------------------------------|-------------------------|--------------------------------------|
| 1. Power indicator                                    | 2. Power Switch         | 3. Brush rotary direction select key |
| 4. Start working indicator                            | 5. Start/Pause          | 6. Pause indicator                   |
| 7. Clockwise rotary working status indicator light    |                         |                                      |
| 8. Inverse hour rotary working status indicator light |                         |                                      |
| 10. Energy decrease key                               | 11. Energy increase key | 12. Time increase key                |
| 13. Time decrease key                                 | 14. Time display        |                                      |


## ●Introduction to Back Panel :

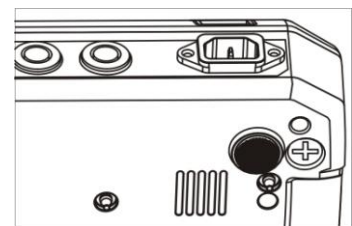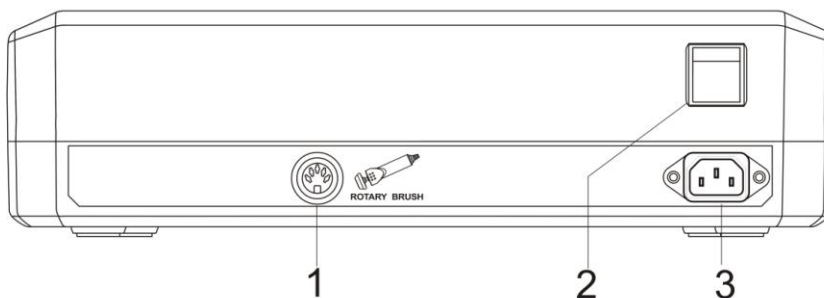

- |                         |                         |
|-------------------------|-------------------------|
| 1.Brush function output | 2. Main power switch    |
| 3. Power input plug     | 4. Protective tube plug |

## ●Connection Instruction :

1. Connect the power cable with host machine;  
The flexible cables of brush handle should be corresponding with outputs at back of host machine.
2. Pick up brushes and plug them to the brush handle (Brushes

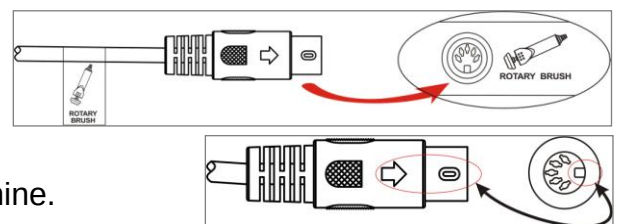

are totally 5 pieces, it can be used in different position ) .

3. Switch on the power at back of host machine, the indicator of main power is lighting.

Press the key " " on control panel,  indicator lights, host machine is stand-by.

Stand-by:  Function indicator lights;

Time display "15";

Energy display "5";

Clockwise rotary working status indicator  lights.

" " Pause indicator at right side  lights.

4. Time setting: press key " " and " " for "TIMER" to set up the working time of host machine, time range 1~30 min.

Energy adjustment: press " " and " " for "ENERGY" to control the energy output of host machine, adjustment levels 1~9.

Select brush rotary direction: Press  to select rotary output status of brush. 

Lights means clockwise rotary output status; and  lights means inverse hour rotary output status; If it is needed to change the rotary direction at the working time, the host machine will come to protection status, and it will transfer automatically after the brush stop to rotaries.

Start or pause output: under the status of stand-by or adjustment completes, press " " to start or pause the output from host machine; when the host is starting output, the start up working indicator  lights; when the host is pause or stand-by, the pause indicator  lights.

5. The prompt volume can be adjusted according to the using environment.

Adjustment: under status of stand-by, press " " key to turn off all displays on control panel (under condition of all indicators no light), press " " key for 3 seconds to enter the function of volume adjustment; Energy display "5" means default volume; press "Energy increase key"  to turn up the voice, and "Energy decrease key"  to turn down the volume,

there are total 9 levels for volume adjustment; no any operation after you have adjusted the suitable volume, the host machine will automatically save the current setting after energy display flashes for 5 seconds, and return to the current status.

#### **Change brush:**

1. Unload brush: Hold the brush handle in one hand, and hold the plastic part of the brush by another hand, then pull in opposite direction, and then unload the brush easily by hearing “click click” sound.
2. Setting up brush: Hold the brush handle in one hand, and hold the plastic part of the brush by another hand, then rotate and push in opposite direction, and then set up the brush easily by hearing “Ka” sound.
3. For unload or set up brush, only can be push or pull the plastic part by hand, but not to do this to the material of brush (e.g. EVA or fuzz etc.). Otherwise will make bad effect to work, even destroy brush.

#### **Cleaning skin:**

1. Firstly use sprayer or clean water softens skin.
2. Wipe skin care production (e.g. cleaning milk or cleaning jelly etc.)equably on skin, according to the characteristic and need to cleaning skin area to choose suitable brush and set up on the brush handle, dip a little water on brush and move it to needed to cleaning skin, turn on host machine then can be come to cleaning and massage.
3. Keep skin moist during the whole cleaning operation.

#### **●Using attention and notices of Brush function:**

1. Before using this brush to work, please check if brush was insert in the brush handle firstly.
2. Please turn down the rotary speed of brush before touching customers’ skin; when touch customers’ skin, it must be slightly to avoid causing customer feel uncomfortable.
3. When use brush to cleaning or massaging skin, please use thickness skin care products and turn down rotary speed of brush, in case of skin care products splashed by fast rotary speed and causing waste and pollution.
3. If it is needed to change brush during skin caring, please stop outputting of the host machine. And wait brush stop to unload and setting up.
4. Please clean and sterilize all the used brushes after finishing work, and keep them in the dry place.

5. It cannot clean brush handle by water but wipe it clean by towel and keep it in the dry place.

### ●Daily Maintenance:

1. Please clean and sterilize all the used brushes after finish working.
2. The instrument should be wiped with dry dishcloth, DO NOT wet in water.
3. Please switch on the main power after using of instrument every day, and takes away the plug when finishes a working day.
4. The instrument should NOT be placed in sunlight, raining, moisture, dust and smoke accumulated area.

### ●Warning :

#### When using Rotary brush function:

1. Do not use it for a very long time, otherwise it will cause brush handle very hot.
2. Do not change brush during it is working.
3. Do not make brush
4. Do not use crystallized brush which has abrasion.
5. During working, the cables between accessory and host machine shall not be drawn excessively.
6. Broken mains cord shall be replaced with specific cord or specific assembly from its manufacturer or maintenance departments. It is strongly recommended that the accessory from original factory shall be used.
7. The host machine shall be turned off and the attaching plug shall be pulled out if any abnormal condition occurs. You shall ask the agency or specific serviceman to check and maintain. The main case is not allowed to be opened.
8. It is forbidden to place or use it at bathroom of wet places.
9. It is forbidden to put it together with inflammable and explosive objects.


### ●Check list for common troubles:

| S/N | Troubles                                                                                                                             | Causes                                    | Remedies                                                                                                                  | Settlements                                                                                  |
|-----|--------------------------------------------------------------------------------------------------------------------------------------|-------------------------------------------|---------------------------------------------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------|
| 1   | When the general supply at the back of the host machine is turned on, the power indicator is off and pressing "POWER" does not work. | The host machine has not been electrified | 1. To check the protective tube of power supply.<br>2. To check the flexible tube of power supply is well plugged or not. | 1. To change the protective tube if necessary.<br>2. To plug the cord of power supply again. |

|   |                                                                                                                                                           |                                                                                                                                       |                                                                                                                                                                                                                     |                                                                                                                                                                                                                                                                                         |
|---|-----------------------------------------------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| 2 | When the general supply is turned on, the light works, but when the host machine is moved, the light glints and the machine cannot be turned on normally. | Power input does not touch well                                                                                                       | <ol style="list-style-type: none"> <li>1. To check the plug of the flexible tube of power supply to see if it is connected well.</li> <li>2. Whether the flexible tube of power supply is broken or not.</li> </ol> | <ol style="list-style-type: none"> <li>1. Poor connection: to plug the connecting plug between the host machine and electric outlet again.</li> <li>2. If the flexible tube of power supply is broken, related parts shall be changed in specific maintenance center nearby.</li> </ol> |
| 5 | The host machine starts output, brush handle does not work.                                                                                               | <ol style="list-style-type: none"> <li>1. The cable of brush does not plug well.</li> <li>2. The cable of brush is broken.</li> </ol> | <ol style="list-style-type: none"> <li>1. To check the cable plug if it is connected well.</li> <li>2. Whether the cable of brush handle is broken or not.</li> </ol>                                               | <ol style="list-style-type: none"> <li>1. Re-plug the cable plug on the back of the host machine.</li> <li>2. Go to the appointed specific maintenance center to change broken brush cable.</li> </ol>                                                                                  |
| 6 | Rotary brush loose                                                                                                                                        | Brush does not set up well                                                                                                            | To check whether the brush was inserted well in the brush handle.                                                                                                                                                   | Take the plug out and re-plug it tight; If re-plug does not work, go to the appointed specific maintenance center to change related accessories.                                                                                                                                        |

### To change protective tube:

Disconnect the power of the host machine first. The cover of protective tube shall be unscrewed by a screwdriver. Check the protective tube inside the protective tube case near to the power supply socket of the host machine, and change it with a new 250V/4A protective tube if it is damaged. (As shown in Fig.: 1. Electric outlet

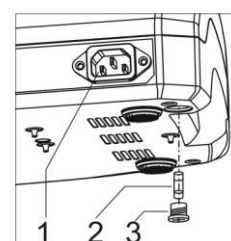

2. Protective tube

3. The cover of protective tube).