

Spot Removal Beauty Instrument User Manual 124222



Dear customers,

Thank you for choosing our product! Please go through this manual before use in order to ensure the maximum performance of the instrument and the safety in use. Please keep properly for further requires.

Function


Microcomputer spot removal machine is the high-tech beauty equipment, with the help of the high frequency current generated by the machine; you can remove freckles, age spots, polyps, etc instantly without bleeding, pain and any side effects.

Specification

1. Easy and fast operation by adopting latest microcomputer technology
2. Smart and reliable treatment by adopting 10 adjustable levels
3. Various methods available by adopting pulse wave

Operation Panel Description

- | | |
|-------------------|----------------------|
| ① Time Decrease | ② Energy Decrease |
| ③ ON/OFF | ④ Time Increase |
| ⑤ Energy Increase | ⑥ Stand-by Indicator |
| ⑦ Time Display | ⑧ On Indicator |

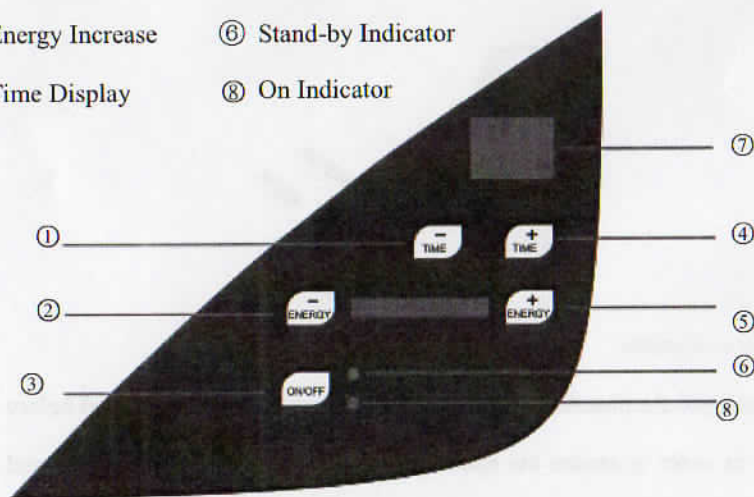

Operation

1. Insert the pin into the handle and insert the handle into the interface.
2. Plug in, then you can hear the startup sound "beep", ⑦ is displayed as 30 (minutes); energy is zero; ⑥ is green, which means standby.
3. Press button ③, indicator ⑧ turns green, the second point in ⑦ flashes, and the machine is started.
4. Increase/decrease the energy at 10 levels by slowly pressing button ⑤ / ②. The energy should be not too high to avoid any injury or not too low to weaken the treatment result. Test the function on some metal part before on client.
5. The pin should be moved forth and back on the treatment position.

Beautician should keep the handle balance not shake to avoid burning.

6. After treatment, shut down the power and clean the pin.

※ For Mole Removal: some moles on the surface of skin but some in the deep. Don't cut through and break the surface one time, should burn and remove in small frame way. To deep mole, should spread proper cell cream on it 2 or 3 days before the wounded cures, to avoid inflammation or skin hollow. Do not spread the cream if the wounded inflame.

※ For Caruncle Removal: gently touch the end with pin and kill the caruncle cell, later about 6-7 days, the caruncle drops. If the caruncle is too big, should remove it several times and not finish it in one time to avoid the wounded covers bad.

Notes

1. If beautician does the treatment with sterilized gloves, hands should first be cleaned and disinfected with alcohol
2. The pin should be sterilized before use, clean the treat area with alcohol then spread it with anesthetic 3 minutes

Caution

1. To move this instrument, it is necessary to lift it up and put it down gently; in case of the failure, please do not knock hard on it so as to prevent malfunction caused by the impact.
2. To move this instrument, it is necessary to keep it horizontal so as to avoid vibration and impact.
3. Do not squeeze the instrument or place heavy stuff onto the instrument so as to prevent shell deformation, damage or malfunction.
4. It is necessary to select the corresponding working power supply in accordance with the marks on the instrument body.
5. It is not allowed to use this instrument in the bathroom.
6. To prevent the lightning accident, when it thunders, please shut down the power switch and pull out the power plug.
7. It is not allowed to inject any fluid into the instrument body so as to prevent the instrument from being damaged due to the short circuit.
8. If the instrument will not be used for a long time, it is necessary to pull out the power plug from the socket.

Technical parameters

1. Operating voltage: As label on the machine
2. The power: As label on the machine
3. Packing size: 27*10*27 (L*W*D) cm